




**Anti-GFP Tag Antibody (Mouse Monoclonal)**  
**Cat#G096**  
**Performance Data**

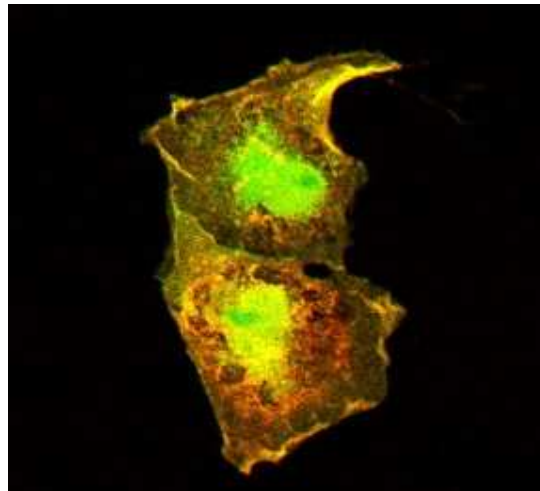

Paraformaldehyde fixed COS-7 cells expressing Myr-N15-PAK2-EGFP construct (Vilas et al.(2006) PNAS 103, 6542). Myr-N15-PAK2-EGFP fluorescence is shown in green. Indirect immunofluorescent detection of N15-PAK-EGFP using G096 monoclonal LGB-1 anti-GFP at 0.05 ug/ml with chicken anti-mouse secondary antibody conjugated to Alexa594 diluted 1/500 is shown in red. Myr-N15-PAK2-EGFP is localized to membrane ruffles and perinuclear vesicular structures (likely Golgi,TGN or late endosomes).